


Titanium Cannulated Screws

A comprehensive compression screw system.

A complete titanium cannulated compression screw system with self drilling/self tapping headed and headless options. Go-EZ™ Headed Cannulated Screws contain Ø 2.0, 2.5, 3.0, 3.5, 4.0, 4.5, 6.5mm and are available in lengths ranging from 6mm to 120mm.† HBS® Headless Cannulated Screws contain Ø 2.5, 3.0mm and are available in lengths ranging from 10mm to 40mm.

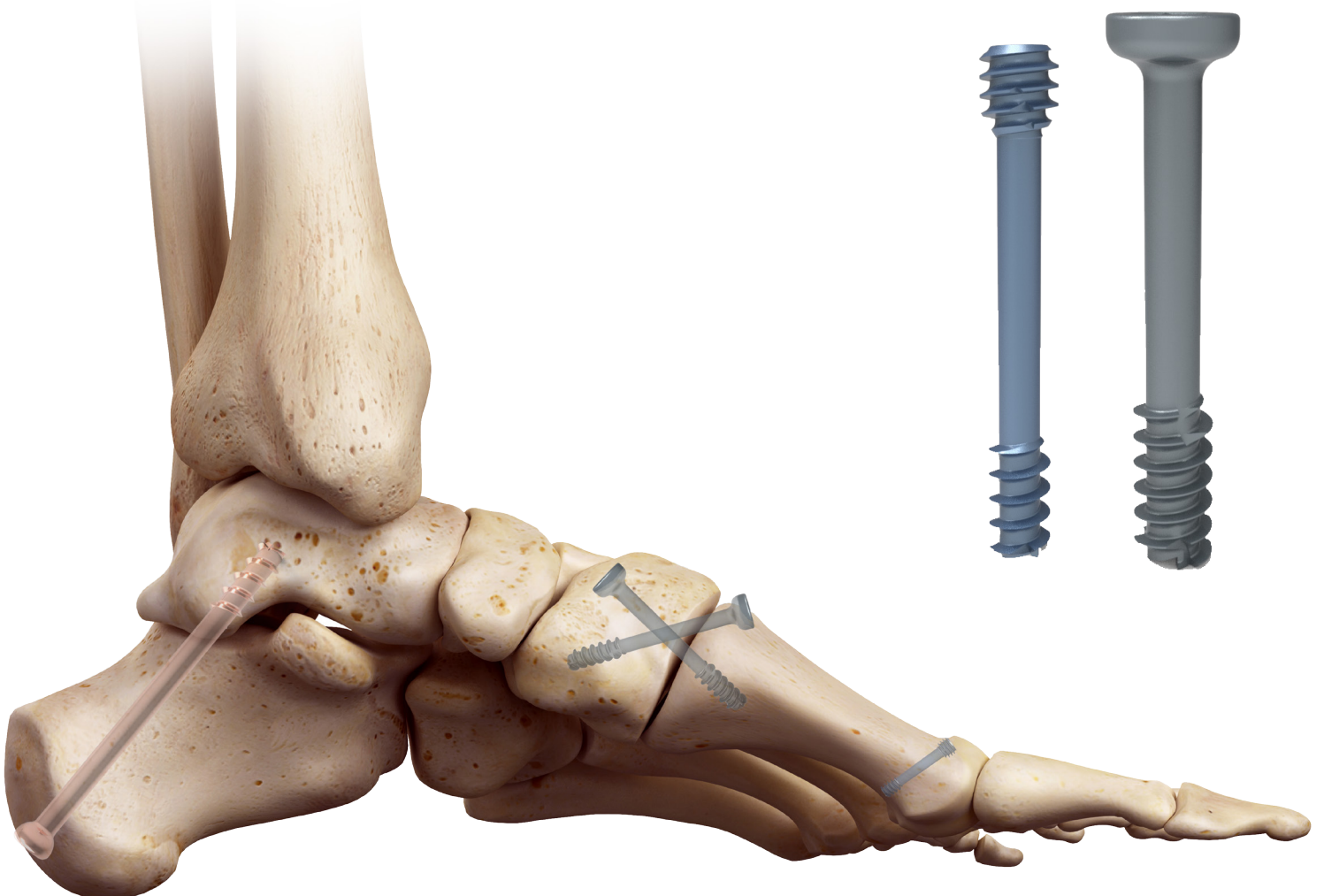

Titanium Cannulated Screws

FINANCIAL FLEXIBILITY

BioPro offers consignment screw inventory and instrument sets. Additionally, our sterile packed screws eliminate the need for bioburden testing and reduce instrument processing costs.

PRODUCT FLEXIBILITY

Screw sizes and quantities are customizable to the surgeon's and/or hospital's needs or preferences.

PROCEDURE FLEXIBILITY

Headed and Headless designs in multiple diameters, lengths and features allow a surgeon to cover a wide variety of procedure options.

Visit www.bioproimplants.com or call us at (810) 982-7777

Go-EZ™ Headed Screw System

Ø 2.0mm | 2.5mm | 3.0mm | 3.5mm | 4.0mm | 4.5mm | 6.5mm



Go-EZ screws have multiple cannulation diameters and feature a hex receiver.




The Go-EZ system features three low profile head designs. The Ø2.0/2.5mm have a washer style head, the Ø3.0,3.5,4.0,4.5mm screws have a standard head design and the Ø6.5mm incorporates a rounded head design.

HBS® Headless Screw System

Ø 2.5mm | 3.0mm



HBS screws have an .035" cannulation diameter and feature a star receiver.

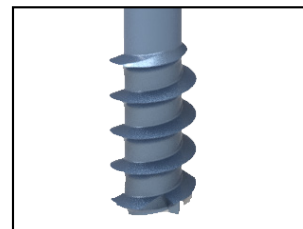

The HBS and Go-EZ screws feature self drilling/self tapping cutting teeth that provide for easy insertion.

† Lengths vary between sizes.



+1-810-982-7777
 info@bioproimplants.com
 2929 Lapeer Road
 Port Huron, MI 48060

