


Foot and Ankle Solutions

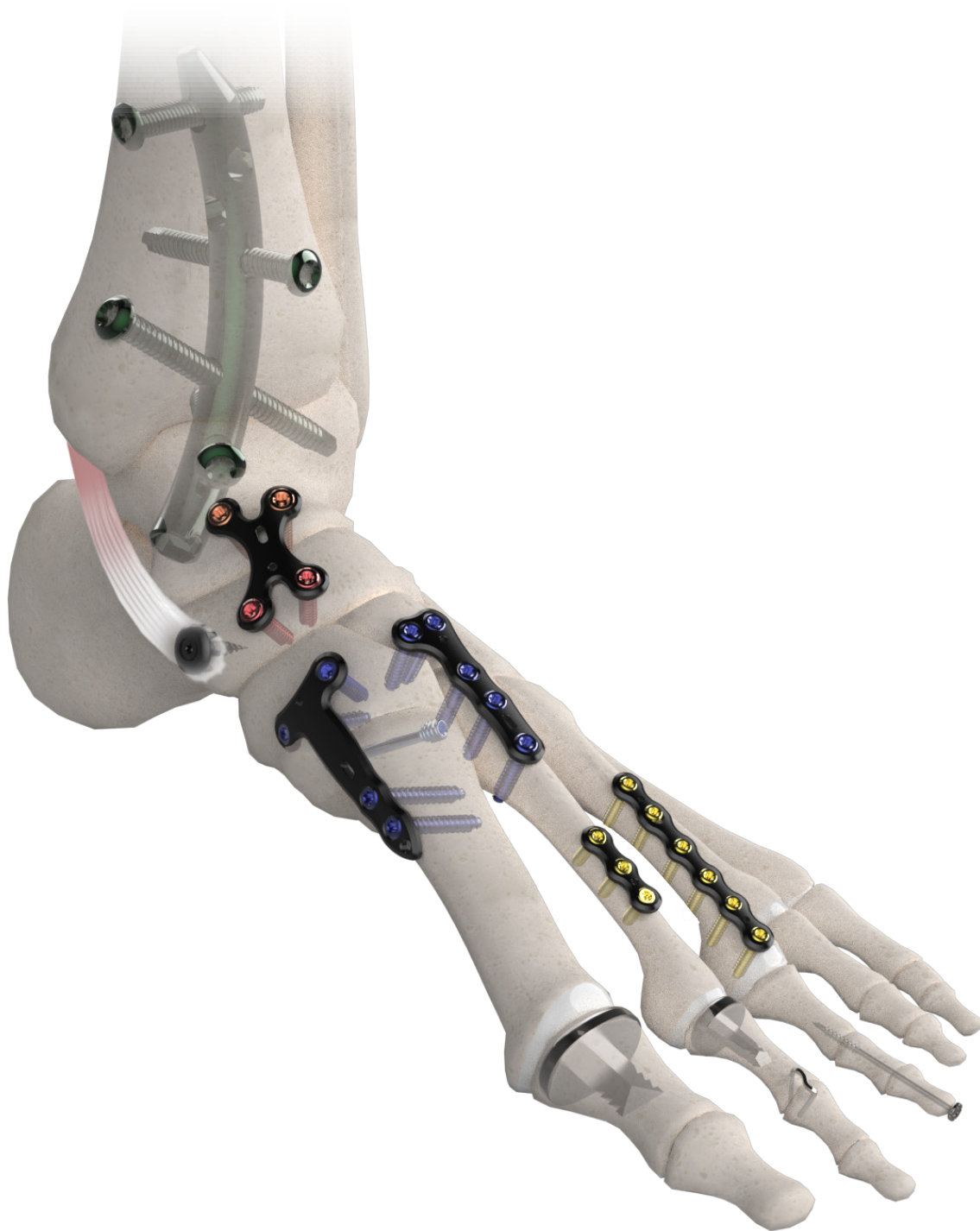

Contact information

Website:
www.bioproimplants.com

Email:
info@bioproimplants.com
orders@bioproimplants.com

Phone:
(810) 982-7777

Fax:
(810) 982-7794

Address:
2929 Lapeer Rd,
Port Huron, MI 48060

Procedure Overview

Hallux Valgus

Metatarsal Osteotomy

Accu-Cut Osteotomy Guides	1
Cannulated Compression Screws	10-11
VIP Tool Screws	12

Akin Osteotomy

Cannulated Compression Screws	9-10
Memory Staple	8
VIP Tool Screws	12

Basal Opening Wedge Osteotomy

Foot Plating System	13
---------------------	----

Lapidus

Cannulated Compression Screws	10-11
Foot Plating System	11-12
Memory Staple	8

Flatfoot

Subtalar Arthrodesis

Cannulated Compression Screws	10-11
-------------------------------	-------

Subtalar Arthroereisis

Horizon Subtalar Implant	6
--------------------------	---

Lateral Column Lengthening (Evans)

Foot Plating System	13
OsteoSinter Wedges	15

Cotton Procedure

OsteoSinter Wedges	15
--------------------	----

Medial Calcaneal Slide

Foot Plating System	13
Cannulated Compression Screws	9-10

Hallux Rigidus

MTPJ Hemiarthroplasty

First MPJ Hemi Implant	2
Silk Toe	3
HemiEDGE Implant	4

MTPJ Arthrodesis

Cannulated Compression Screws	10-11
Foot Plating System	13
Memory Staple	8

Lesser Toe Deformities

DIPJ/PIPJ Arthrodesis

Cannulated Compression Screws	10-11
Digital Compression Screws	8

MTPJ Hemiarthroplasty

Lesser MPJ Hemi Implant	5
-------------------------	---

Fractures

Metatarsal Fractures

Cannulated Compression Screws	10-11
Foot Plating System	13

Tendon Reattachment

Tendon Anchor System	6
----------------------	---

Ankle Arthritis

Shotel Ankle Arthrodesis Nail	14
-------------------------------	----

Accu-Cut® Osteotomy Guide System

The Accu-Cut® Osteotomy Guide System is a disposable, sterile packaged system that includes guide(s), two universal saw blades, and two 0.045" double trocar K-wires. Guides are available for the most popular distal metatarsal osteotomies to address mild to moderate hallux valgus deformity.

Using an Accu-Cut® Osteotomy Guide ensures precise and accurate cuts every time which is proven to provide improved stability and reduce the risk of displacement.

- All-inclusive sterile package provides the surgeon with everything needed to perform the desired osteotomy.
- Consistent and accurate correction reduces the risk of displacement.
- Simple technique allows performance of complex osteotomies quickly and efficiently.



Chevron (Austin) Guide

The standard Chevron Guide offers a classic v-shape chevron osteotomy to address hallux valgus by lateral transposition of the metatarsal head.

Youngswick Guide

The Youngswick Guide is available in either 1mm, 2mm, or 3mm dorsal wedges for performing a Youngswick chevron modification to address both hallux valgus and hallux rigidus.



Long-Arm Chevron Guide

Z (Scarf) Guide

DMAA Guide



2-Stage Guide

Item #	Description
19505	Austin/Chevron Guide
19506	Youngswick 1mm
19507	Youngswick 2mm
19508	Youngswick 3mm
19510	Long Arm Chevron Guide
19511	Z (Scarf) Guide
19720	DMAA Guide
19509	2-Stage PASA/DMAA Guide

First MPJ Hemi Implant

The First MPJ Hemi Implant is the only hemi backed by 70+ years of clinical data. Our low-profile, press-fit implants are designed to replace the articular surface of the proximal phalanx in a painful, arthritic metatarsophalangeal (MTP) joint. The procedure offers restoration of motion and pain relief, with a minimal bone resection technique.^{1,2}

- 20+ year survivorship data.¹
- +95% implant survivorship on average.*
- 97% of patients would recommend the procedure.²
- Immediate weight-bearing.
- ≤6 weeks return to activities.²

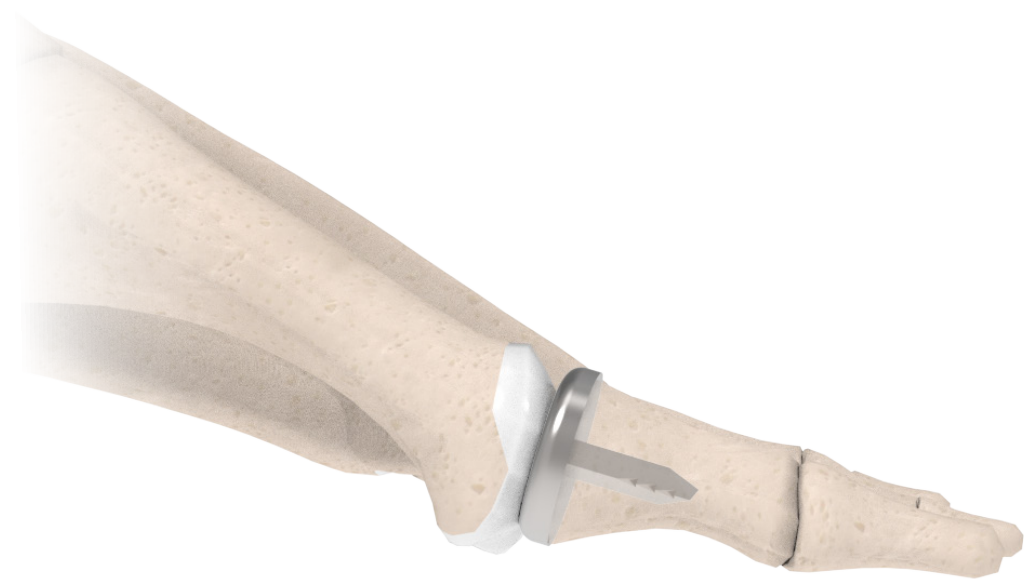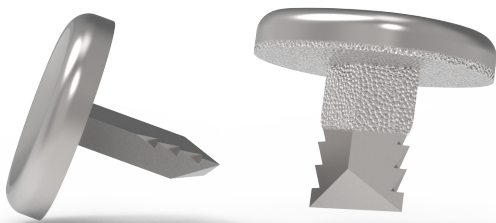


	SM 17mm	M/S 18.5mm	MD 20mm	M/L 21.5mm	LG 23mm
Cobalt Chrome PC	10412	17034	10413	14960	10414
Cobalt Chrome NPC	10060	17033	10061	14958	10062
Titanium PC	17035	17197	17036	17037	17038

REFERENCES
1. Townley, MD, Taranow, DO. A metallic hemiarthroplasty resurfacing prosthesis for the hallux metatarsophalangeal joint. Foot & Ankle International 1994;15(11):575-80
2. Beekhuizen, Stefan R. et al. Long-Term Results of Hemiarthroplasty Compared With Arthrodesis for Osteoarthritis of the First Metatarsophalangeal Joint. The Journal of Foot and Ankle Surgery , Volume 57 , Issue 3 , 445 - 450
* Data on file

SilkToe

The BRM Silktoe® Implant is a new generation flexible silicone spacer for first metatarsophalangeal (1st MTP) joint arthroplasty.


- Available in five sizes in both left and right orientations to preserve valgus alignment
- Anatomically contoured hinge to allow for improved range of motion
- Features an angled resection that preserves the flexor hallucis brevis tendon
- Titanium grommets to allow for ingrowth and improved stability

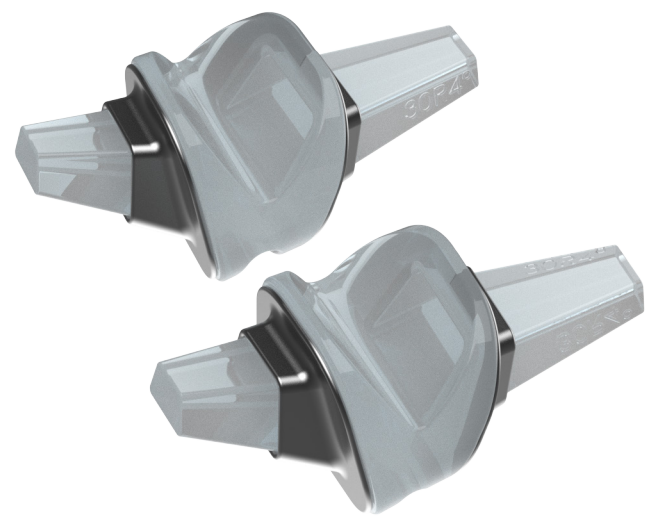


Item #	Description
BRDNEWP2020L	Silktoe Spacer - Size 20 Left
BRDNEWP2020R	Silktoe Spacer - Size 20 Right
BRDNEWP3030L	Silktoe Spacer - Size 30 Left
BRDNEWP3030R	Silktoe Spacer - Size 30 Right
BRDNEWP3535L	Silktoe Spacer - Size 35 Left
BRDNEWP3535R	Silktoe Spacer - Size 35 Right
BRDNEWP4040L	Silktoe Spacer - Size 40 Left
BRDNEWP4040R	Silktoe Spacer - Size 40 Right
BRDNEWP5050L	Silktoe Spacer - Size 50 Left
BRDNEWP5050R	Silktoe Spacer - Size 50 Right

HemiEDGE™

The HemiEDGE is a simple, low-profile, press-fit implant designed to replace the articular surface of the proximal phalanx in a painful, arthritic metatarsophalangeal (MTP) joint. Based on the clinical success of our First MPJ Hemi Implant, the HemiEDGE, incorporates an overlapping edge extending around the medial, lateral and dorsal aspects of the implant. Partially encompassing the cortex of the phalanx helps ensure proper implant sizing, improve implant stability, and reduce the potential of bony overgrowth. The ultra thin profile of 2.5mm requires minimal bone resection at the base of the proximal phalanx for proper implantation.

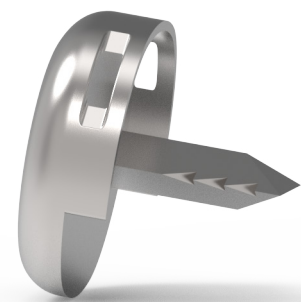

- Provides additional stability and ensures proper sizing.
- Reduces the risk of potential bone over-growth.
- Maintains low 2.5mm profile for minimal bone resection.

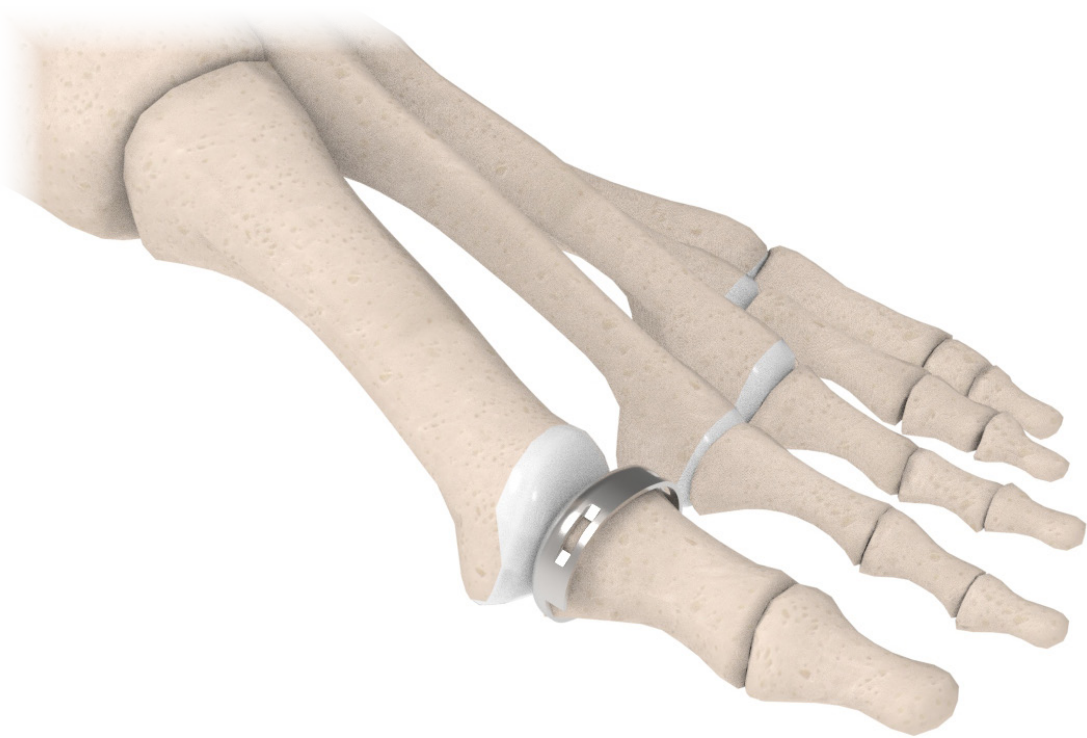

	SM 17mm	M/S 18.5mm	MD 20mm	M/L 21.5mm	LG 23mm
HemiEDGE	19538	19539	19540	19541	19542

Lesser MPJ Implant

The Lesser MPJ Hemi Implant is a low profile, press-fit implant designed to replace the articular surface of the proximal phalanx in a painful, arthritic lesser metatarsophalangeal (MTP) joint. The implant offers restoration of motion and pain relief, with a minimal bone resection technique. The implants are manufactured from cobalt chrome and available in 9 sizes. The implant is cannulated which allows a k-wire to be temporarily passed through the metatarsal head to stabilize the entire digit when performed in conjunction with a hammertoe correction.



- Allows for correction of lesser toe deformities such as Freiberg's infraction and osteoarthritis, with immediate weight-bearing.
- May be used in conjunction with hammertoe correction.
- Low profile implant only requires 2mm of resection on the proximal phalanx.





	8mm	8.75mm	9.5mm	10.25mm	11.00mm	11.75mm	12.50mm	13.25mm	14.00mm
Cannulated	17324	17325	17326	17327	17328	17329	17330	17331	17332

Horizon® Subtalar Implant

The Horizon® Subtalar Implant is indicated for use in the treatment of the hyper-pronated foot and stabilization of the subtalar joint. It is designed to block forward, downward and medial displacement of the talus, thus allowing normal subtalar joint motion but blocking excessive pronation. The implant incorporates a smooth, slightly tapered lateral trailing edge; assuring that the talus makes contact with a smooth, flat surface when the foot is pronated, not a sharp thread.

- Minimally invasive and reversible procedure.
- Allows for quick recovery and pain relief compared to other surgical interventions.
- Smooth lateral surface ensures both the talus and calcaneus interface with a smooth surface during pronation.
- Slight taper stabilizes the implant during pronation and protects against lateral migration.

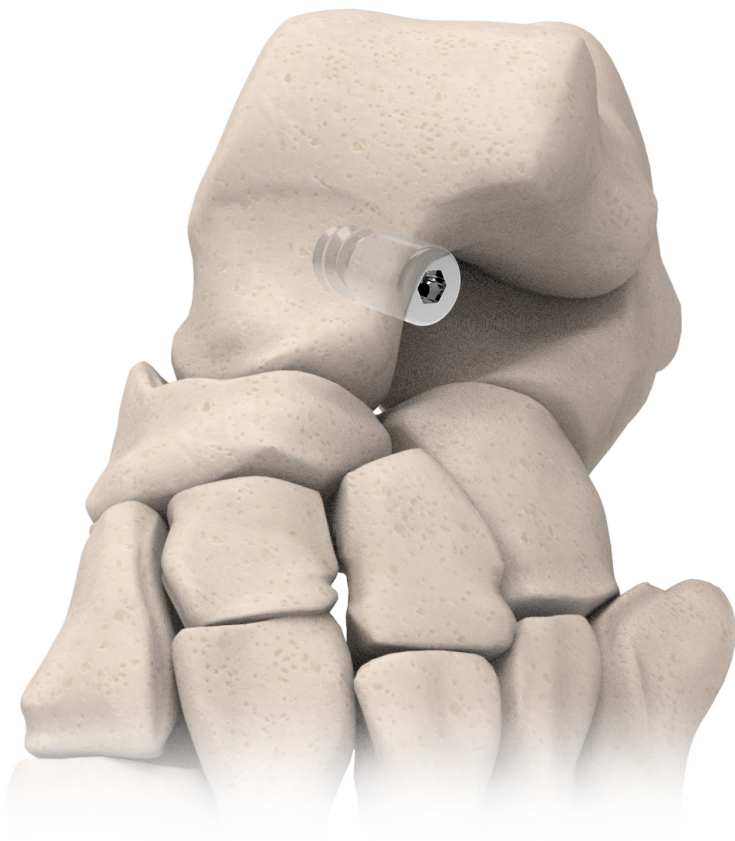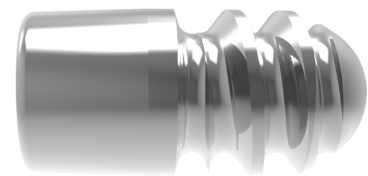



	6mm	8mm	9mm	10mm	11mm	12mm
Standard	17478	17221	17222	17223	17224	17225
Hybrid	N/A	17085	17086	17087	17088	17089

1. Ozan F, Doğar F, Genç K, Koyuncu, Vatansever F, Duygulu F, Altay T. Symptomatic flexible flatfoot in adults: subtalar arthroereisis. Ther Clin Risk Manag. 2015;11:1597-1602 <https://doi.org/10.2147/TCRM.S90649>

Tendon Anchor System™

The Tendon Anchor System (T.A.S.) is a toothed, titanium anchor designed for soft tissue reattachment to bone. The sutureless design simplifies soft tissue reattachment and increases pull-out strength.¹ The extremely low profile design allows the anchor to recess into the tendon preventing implant prominence, while still allowing easy removal if necessary.

- No suture required, eliminating the tendon-suture interface, the most common failure mode in soft tissue attachment.²
- Simple place and screw technique with no pre-drilling required, reducing OR time.[†]
- Provides greater surface area coverage resulting in 3 times greater pull-out strength vs standard suture anchors.¹
- Low-profile (1mm) anchors reduce soft tissue irritation.
- No tendon detachment required for implantation or removal.
- Blunt teeth prevent over-compression allowing for microcirculatory blood flow crucial to the healing process.

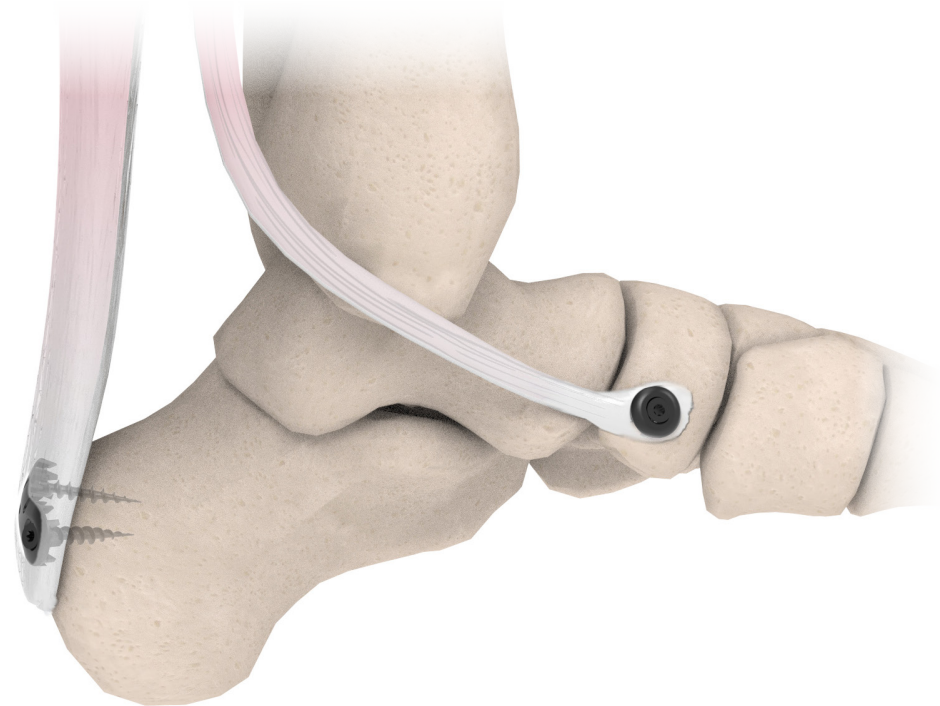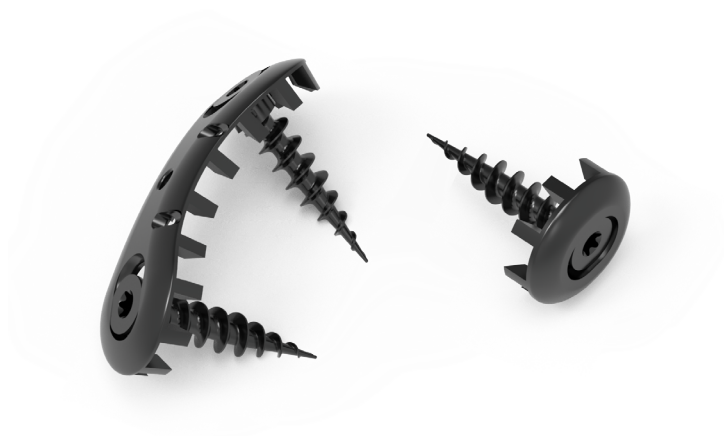


Item #	Description	Size
20070	One-Hole Tendon Anchor System	10MM
20159	Two-Hole Tendon Anchor System	26MM

Memory Staple

The BioPro Memory Staple is a 2-leg nitinol memory-alloy staple designed for bone fixation. The unique S-Bend bridge of the Memory Staple ensures even compression across the fusion site, while maintaining a low profile against the bone.

The Memory Staple is individually sterile packaged and available in 17 standard sizes. The packaging features a color coding system that coordinates staple leg size with drill diameter. The system offers standard legs or offset legs. Standard leg staples are utilized when a flush surface exists while the offset staples offer different leg lengths to accommodate uneven bone surfaces.

- 17 standard sizes with varying bridge widths, leg lengths, and wire sizes, allowing fixation in most hand or foot procedures
- Easy to use instrumentation for fast and direct visual insertion may reduce procedure time
- Body-heat activated design offers delayed compression, ensuring proper seating and positioning of the staple prior to compression

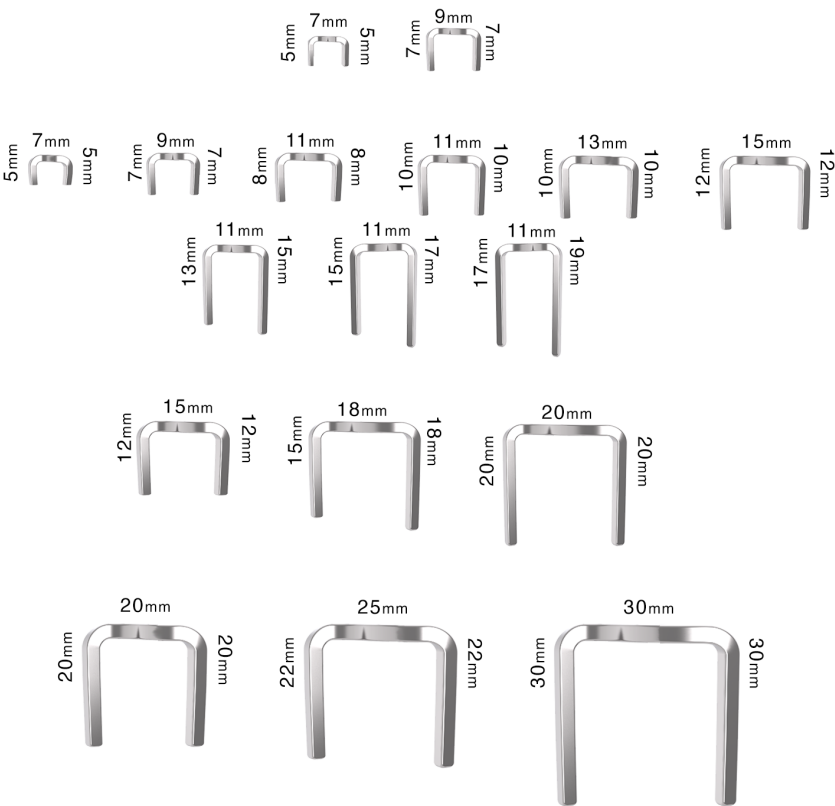

Item #	Bridge Width	Leg Length	Wire Size
17637	7mm	5mm	1.2mm
17628	7mm	5mm	1.5mm
17638	9mm	7mm	1.2mm
17629	9mm	7mm	1.5mm
17630	11mm	8mm	1.5mm
17631	11mm	10mm	1.5mm
17634	11mm	15mm/13mm	1.5mm
17635	11mm	17mm/15mm	1.5mm
17636	11mm	19mm/17mm	1.5mm
17632	13mm	10mm	1.5mm
17633	15mm	12mm	1.5mm
17625	15mm	12mm	2.0mm
17626	18mm	18mm/15mm	2.0mm
17627	20mm	20mm	2.0mm
17622	20mm	20mm	2.0x3.0mm
17623	25mm	22mm	2.0x3.0mm
17624	30mm	30mm	2.0x3.0mm

Digital Compression Screw

The Digital Compression Screw (DCS) is a solid, stainless steel screw specifically designed to address digital fusions. The 1.5mm and 1.8mm diameters range from 20mm to 55mm in length, allowing for fusion of the DIPJ, PIPJ, or both. The screw features a 1.1mm shaft diameter with either 1.5mm or 1.8mm threads for a lag screw design. Our unique over-drilling technique allows compression and easy removal post fusion.



- Lag design combined with a unique over-drilling technique allows for compression across the fusion site
- Threads purchase into cortical bone, improving stability and reducing the chance of rotation or residual flexion
- Unlike intramedullary fixation devices, the DCS allows for fusion of both the DIPJ and PIPJ reducing the chance of mallet toe
- Buried beneath the skin, it allows patients to wear a stiff-soled postoperative shoe
- Designed for simple in-office screw removal post fusion



Item #	Description	Length
17201	Digital Compression Screw 1.5mm	20mm
17202	Digital Compression Screw 1.5mm	25mm
17110	Digital Compression Screw 1.5mm	30mm
17111	Digital Compression Screw 1.5mm	35mm
17112	Digital Compression Screw 1.5mm	40mm
17113	Digital Compression Screw 1.5mm	45mm
17114	Digital Compression Screw 1.5mm	50mm
17115	Digital Compression Screw 1.5mm	55mm
17203	Digital Compression Screw 1.8mm	20mm
17204	Digital Compression Screw 1.8mm	25mm
17116	Digital Compression Screw 1.8mm	30mm
17117	Digital Compression Screw 1.8mm	35mm
17118	Digital Compression Screw 1.8mm	40mm
17119	Digital Compression Screw 1.8mm	45mm
17120	Digital Compression Screw 1.8mm	50mm
17121	Digital Compression Screw 1.8mm	55mm

Go-EZ™ Headed Compression Screws

The Go-EZ™ Headed Cannulated Compression Screw System features screws from 2.0mm to 6.5mm in diameter. The screws are manufactured from implant grade Titanium and are individually sterile packaged.

The system features low profile heads minimizing implant prominence and self-drilling and self-cutting tips for effortless insertion without pre-drilling.

Our screw sizes and quantities are customizable to the surgeon's/hospital's needs or preferences.



Go-EZ Cannulated Screws								
Length (mm)	Ø2.0 (Blue)	Ø2.5 (Purple)	Ø3.0 (Gold)	Ø3.5 (Green)	Ø4.0 (Grey)	Ø4.5 (Brown)	Length (mm)	Ø6.5 (Rose)
6	18162	18172	18185	18203	18226	18239	30	19651
8	18163	18173	18186	18204	18227	18240	32	19652
10	18164	18174	18187	18205	18228	18241	34	19653
12	18165	18175	18188	18206	18229	18242	36	19654
14	18166	18176	18189	18207	18230	18243	38	19655
16	18167	18177	18190	18208	18231	18244	40	19656
18	18168	18178	18191	18209	18232	18245	45	19657
20	18169	18179	18192	18210	18233	18246	50	19658
22	18170	18180	18193	18211	18234	18247	55	19659
24	18171	18181	18194	18212	18235	18248	60	19660
26	19265	18182	18195	18213	18236	18249	65	19661
28	19266	18183	18196	18214	18237	18250	70	19662
30	19267	18184	18197	18215	18238	18251	75	19663
32	19268	19250	18198	18216	19041	19031	80	19664
34	19269	19251	18199	18217	19042	19032	85	19665
36	19270	19252	18200	18218	19043	19033	90	19666
38	19271	19253	18201	18219	19044	19034	95	19667
40	19272	19254	18202	18220	19045	19035	100	19668
42	19273	19255		18221	19046	19036	105	19669
44	19274	19256		18222	19047	19037	110	19670
46	19275	19257		18223	19048	19038	115	19671
48	19276	19258		18224	19049	19039	120	19672
50	19277	19259		18225	19050	19040		

*Shaded area indicates special order only

HBS® Headless Compression Screws

The HBS® Headless Cannulated Compression Screw System features 2.5mm and 3.0mm screws with various cutting tip styles available. The screws are manufactured from implant grade Titanium and are individually sterile packaged.

The 2.5mm and 3.0mm self-tapping (ST) screws are available from 10mm to 30mm in length with 1mm increments to allow for maximal bone purchase.

The 2.5mm and 3.0mm self-drilling and self-tapping (SDST) screws are available from 10mm to 40mm in length with 2mm increments. The self-cutting teeth eliminate the need for pre-drilling.

The 3.0mm Fully Threaded (FT) screw is available from 10 to 40mm in length with 2mm increments. The fully threaded design allows for downsizing without sacrificing thread purchase.



HBS ST SCREWS		
Length (mm)	Ø2.5 (Purple)	Ø3.0 (Grey)
10	19545	19566
11	19546	19567
12	19547	19568
13	19548	19569
14	19549	19570
15	19550	19571
16	19551	19572
17	19552	19573
18	19553	19574
19	19554	19575
20	19555	19576
21	19556	19577
22	19557	19578
23	19558	19579
24	19559	19580
25	19560	19581
26	19561	19582
27	19562	19583
28	19563	19584
29	19564	19585
30	19565	19586

HBS SDST SCREWS		
Length (mm)	Ø2.5 (Rose)	Ø3.0 (Blue)
10	20500	20516
12	20501	20517
14	20502	20518
16	20503	20519
18	20504	20520
20	20505	20521
22	20506	20522
24	20507	20523
26	20508	20524
28	20509	20525
30	20510	20526
32	20511	20527
34	20512	20528
36	20513	20529
38	20514	20530
40	20515	20531

HBS Fully Threaded	
Length (mm)	Ø3.0 (Silver)
10	22540
12	22541
14	22542
16	22543
18	22544
20	22545
22	22546
24	22547
26	22548
28	22549
30	22550
32	22551
34	22552
36	22553
38	22554
40	22555

*Shaded area indicates special order only

Tool VIP Screw

The TOOL VIP Screw System offers three different diameter screws and specialized instrumentation designed for minimally invasive surgery. The neutral pitch allows surgeons to maintain length of the osteotomy site and optimize stability.

- The screws are manufactured from Titanium Alloy and feature a 45 degree beveled head that sits flush with the bone upon insertion.
- The screws are fully threaded with a neutral pitch and feature self-drilling + self-tapping cutting tips.
- The three different screw diameters allow for use in multiple applications. The screws are color coded and sterile packaged for improved traceability.



Length	2.5MM	3.0MM	4.0MM
16MM	SSL125016S	SSL130016S	SSL140016S
18MM	SSL125018S	SSL130018S	SSL140018S
20MM	SSL125020S	SSL130020S	SSL140020S
22MM	SSL125022S	SSL130022S	SSL140022S
24MM	SSL125024S	SSL130024S	SSL140024S
26MM	SSL125026S	SSL130026S	SSL140026S
28MM	SSL125028S	SSL130028S	SSL140028S
30MM	SSL125030S	SSL130030S	SSL140030S
32MM	SSL125032S	SSL130032S	SSL140032S
34MM	SSL125034S	SSL130034S	SSL140034S
36MM	SSL125036S	SSL130036S	SSL140036S
38MM	SSL125038S	SSL130038S	SSL140038S
40MM	SSL125040S	SSL130040S	SSL140040S
42MM		SSL130042S	SSL140042S
44MM		SSL130044S	SSL140044S
46MM		SSL130046S	SSL140046S
48MM		SSL130048S	SSL140048S
50MM		SSL130050S	SSL140050S
52MM		SSL130052S	SSL140052S
54MM		SSL130054S	SSL140054S
56MM		SSL130056S	SSL140056S
58MM		SSL130058S	SSL140058S
60MM		SSL130060S	SSL140060S



Tool VIP Screws are manufactured by BRM Extremities. For more information and indications for use visit www.brm-extremities.com

Comprehensive Foot Plating System

A comprehensive forefoot, midfoot and hindfoot plating system merges modern technology with plating basics. This system is truly balanced by a comprehensive selection, providing surgeons with plates of all sizes and shapes to match the anatomical structures of the foot while not overloading the tray with unnecessary options.

The 45 universal plates are categorized by screw diameter including 2.0mm plates, 2.8mm plates and 3.0-3.5mm plates. The indication specific plating options address a variety of procedures including MTP arthrodesis, TMT-1 arthrodesis, basal (opening wedge) osteotomy, calcaneal sliding osteotomy, lateral column lengthening (Evans osteotomy), calcaneal fractures and medial column arthrodesis.

Each plate features MVA (multiple variable angle) locking technology that allows screws to lock into the plate up to 25 degrees. The plates range in thickness from 1.3 to 1.8mm, depending on the procedure needs, and feature chamfered edges to minimize soft tissue irritation. Manufactured from grade 4 titanium with type II anodization, the plates have greater wear resistance, higher fatigue strength and improved biocompatibility.*

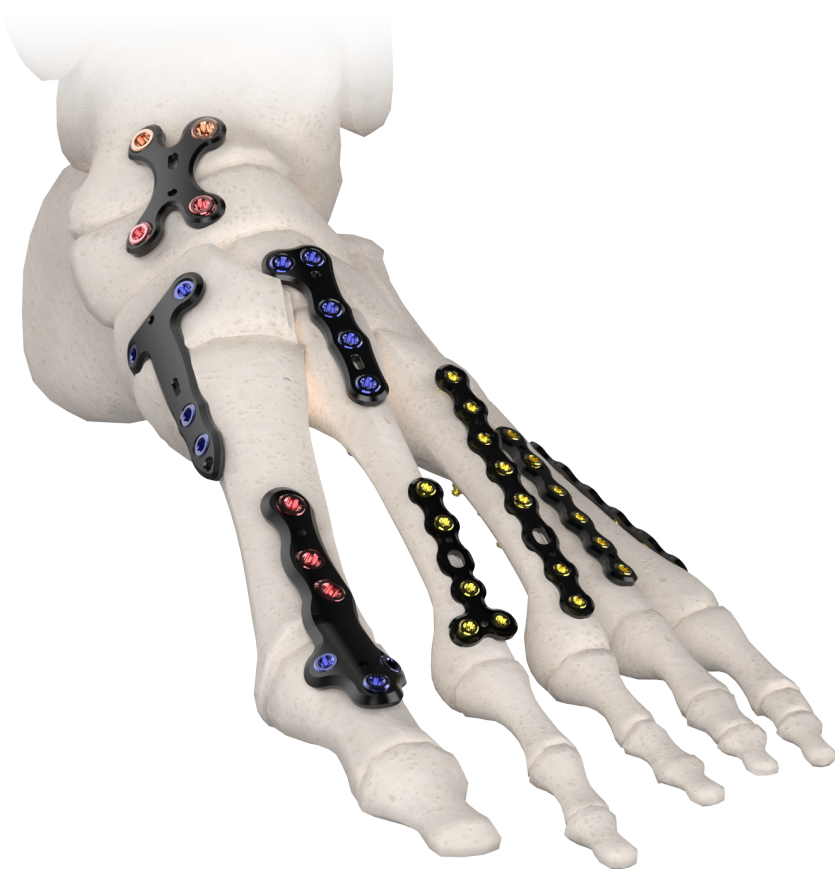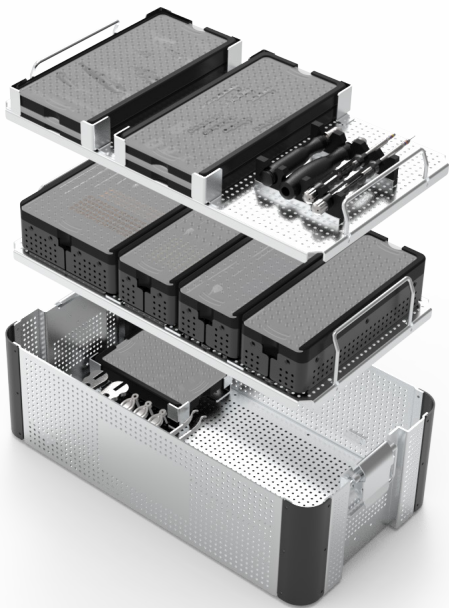


Item #	Description
22929	Base Tray
22926	2.0 Screw + Plate Module
22927	2.8 Screw Module
22930	2.8 Plate Module
22928	3.0/3.5 Screw Module
22931	3.0/3.5 Plate Module
22929	MTP Plate Module
22932	Medial Column Plate Module
22934	Osteotomy Plate Module
22933	Calcaneal Plate Module

*when compared to plates with type III anodization

Shotel™ Ankle Arthrodesis Nail System

The Shotel™ Ankle Arthrodesis Nail System is an intramedullary (IM) nail intended to provide compression and rigid fixation for primary ankle fusions. The unique curved design allows an approach through the medial side of the talus and is designed to achieve fusion at the tibiotalar joint while allowing unrestricted motion to remain at all other joints. The Shotel™ Ankle Arthrodesis Nail System provides a minimally invasive surgical option that may allow early weight bearing.

- The solid titanium IM fusion nail is available in a 9mm diameter in two configurations; left and right. The fusion nail can accommodate up to six bone screws.
- The 5mm self-tapping titanium bone screws are fully threaded and can be used to provide static or dynamic cross-locking through the provided nail holes. The bone screws are available in various lengths from 20–80mm in 2mm increments.

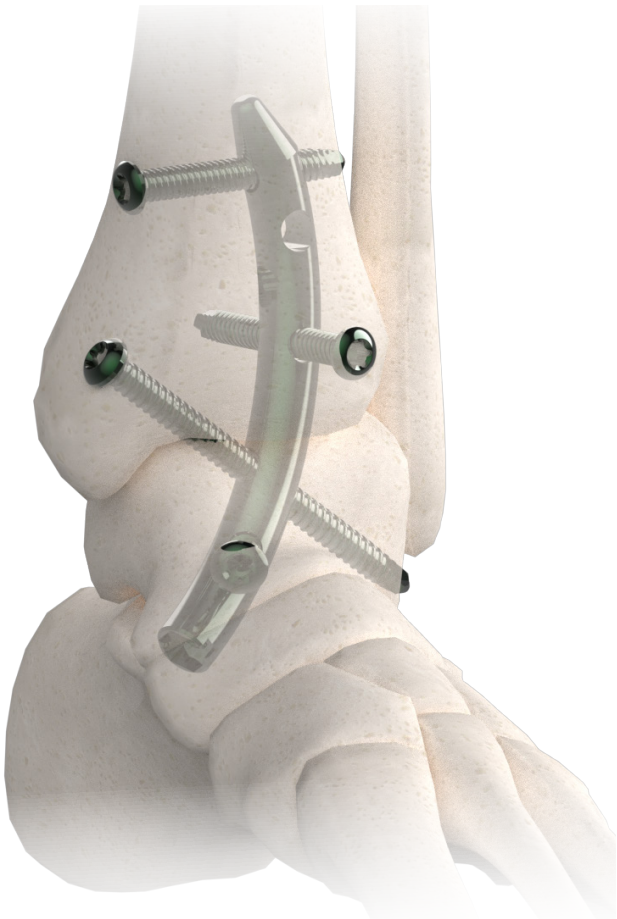

ITEM #	DESCRIPTION
22909	SHOTEL INSTRUMENT KIT COMPLETE
21674	STEINMANN PIN SINGLE TROCAR ROUND 12" 2.4 MM
22687	CURVED GUIDE WIRE SS
22390	FLEXIBLE REAMER WITH MULTIPLE FLUTE TIP, 10.5 MM
21659	4.2 MM DRILL BIT

OsteoSinter® Evans and Cotton Wedges

OsteoSinter® EVANS and COTTON Wedges are porous titanium implants engineered for use in the surgical correction of adult-acquired flatfoot deformities. Designed specifically for Evans and Cotton osteotomies, these implants support structural realignment and biomechanical restoration of the foot.

- Manufactured using a specialized process that replicates the architecture and porosity of cancellous bone, the wedges promote enhanced osseointegration and long-term stability.
- The system includes a comprehensive range of 18 implant sizes—12 for Evans and 6 for Cotton osteotomies—to facilitate precise anatomical correction.
- All implants and instruments are sterile-packaged and ready for immediate, intraoperative use.



Item #	Description
P00209	OsteoSinter® EVANS wedge 18W x 18L x 6H
P00200	OsteoSinter® EVANS wedge 18W x 18L x 8H
P00201	OsteoSinter® EVANS wedge 18W x 18L x 10H
P00202	OsteoSinter® EVANS wedge 18W x 18L x 12H
P00210	OsteoSinter® EVANS wedge 20W x 20L x 6H
P00203	OsteoSinter® EVANS wedge 20W x 20L x 8H
P00204	OsteoSinter® EVANS wedge 20W x 20L x 10H
P00205	OsteoSinter® EVANS wedge 20W x 20L x 12H
P00211	OsteoSinter® EVANS wedge 22W x 22L x 6H
P00206	OsteoSinter® EVANS wedge 22W x 22L x 8H
P00207	OsteoSinter® EVANS wedge 22W x 22L x 10H
P00208	OsteoSinter® EVANS wedge 22W x 22L x 12H
P00401	Single-use OsteoSinter® EVANS instrument kit
P00300	OsteoSinter® COTTON wedge 15L x 4.5H
P00301	OsteoSinter® COTTON wedge 15L x 5.5H
P00302	OsteoSinter® COTTON wedge 15L x 6.5H
P00303	OsteoSinter® COTTON wedge 20L x 4.5H
P00304	OsteoSinter® COTTON wedge 20L x 5.5H
P00305	OsteoSinter® COTTON wedge 20L x 6.5H
P00413	Single-use OsteoSinter® COTTON instrument kit



Call us at 1-810-982-7777 to schedule a case today.

This content is provided as an educational tool only and is not meant as medical advice in the usage of specific BioPro products. A healthcare professional must use their professional judgment in making any final determinations in product usage and technique. The product's Instructions for Use, should always be reviewed prior to surgery. Postoperative management is patient-specific and dependent on the treating professional's assessment. Individual results will vary and not all patients will experience the same postoperative activity level or outcomes. This information does not constitute medical, legal, or any other type of professional advice and should not be relied upon as such. It is not to be redistributed, duplicated, or disclosed without the express written consent of BioPro, Inc.



BioPro, Inc.
2929 Lapeer Road, Port Huron, MI 48060, USA
info@bioproimplants.com | 1-810-982-7777
www.bioproimplants.com



intertek