


HemiEDGE™

The newest advancement in hemi resurfacing implants.

Based on the clinical success of our First MPJ Hemi Implant, the HemiEDGE™ incorporates an overlapping edge extending around the medial, lateral and dorsal aspects of the implant. Partially encompassing the cortex of the phalanx helps ensure proper implant sizing, improve implant stability, and reduce the potential of bony overgrowth.



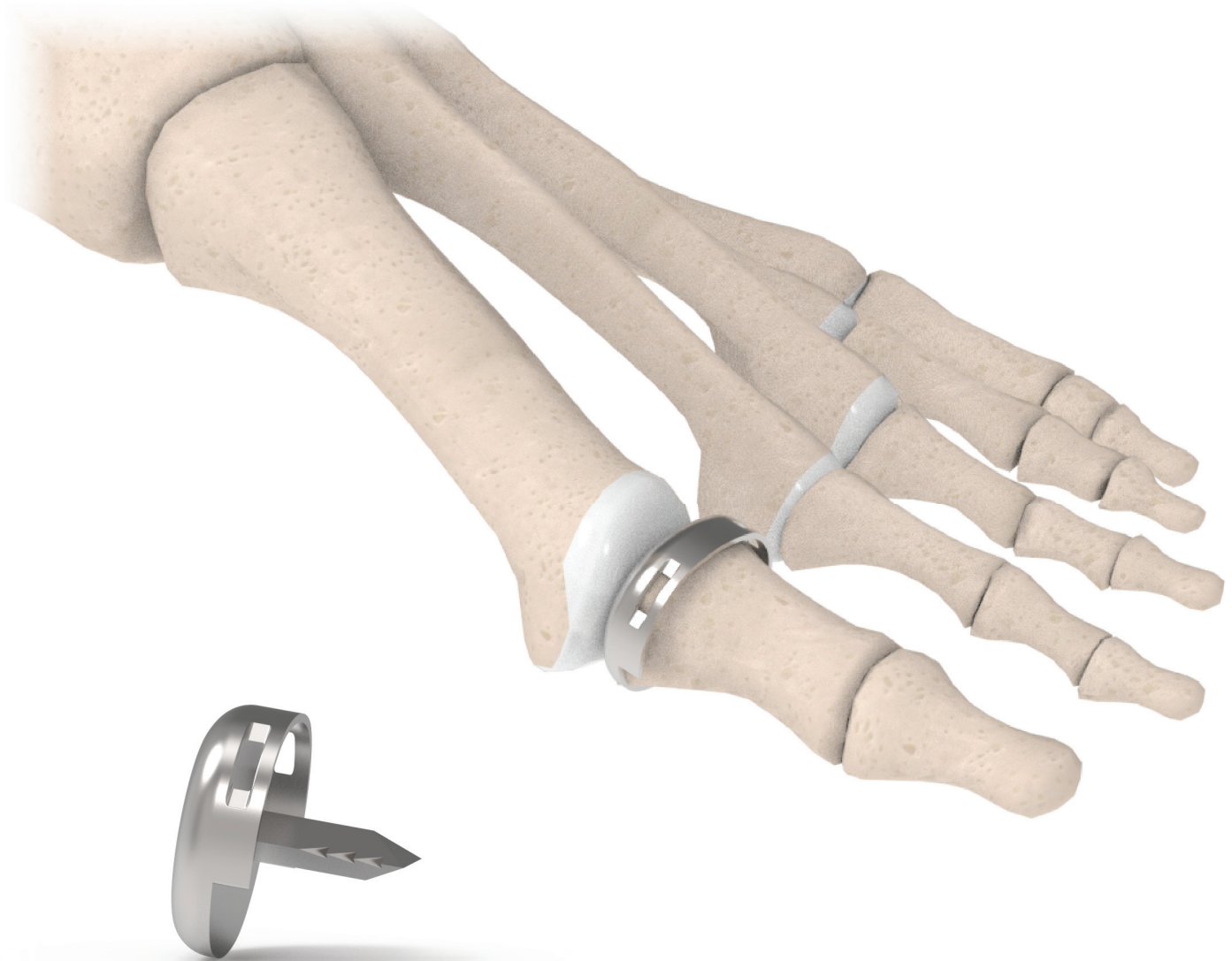

The benefits of the EDGE.

The overlapping Edge is the main design feature differentiating the HemiEDGE from traditional phalangeal implants. This edge extends over the dorsal, medial and lateral cortices beyond the resected base of the phalanx. The plantar surface does not include the edge, preserving the flexor attachment and avoiding any implant prominence on the weight bearing surface. The edge surrounds the phalanx to help support and stabilize the implant. This improves implant stability, assists in proper placement during implantation, and reduces the potential of bony overgrowth.³

Why use BioPro MPJ Hemi Implants?

- ✓ 20+ year survivorship data¹
- ✓ +95% implant survivorship on average*
- ✓ 97% of patients would recommend the procedure²
- ✓ Immediate weight-bearing
- ✓ ≤6 weeks return to activities²

	SM 17MM	M/S 18.5MM	MD 20MM	M/L 21.5MM	LG 23MM
HemiEDGE	19538	19539	19540	19541	19542

REFERENCES

1. Townley, MD, Taranow, DO. A metallic hemiarthroplasty resurfacing prosthesis for the hallux metatarsophalangeal joint. *Foot & Ankle International* 1994;15(11):575-80
2. Beekhuizen, Stefan R. et al. Long-Term Results of Hemiarthroplasty Compared With Arthrodesis for Osteoarthritis of the First Metatarsophalangeal Joint. *The Journal of Foot and Ankle Surgery* , Volume 57 , Issue 3 , 445 - 450
3. Jerome A. Slavitt, DPM, FACFAS . A Closer Look At A First MPJ Hemi-Implant For The Treatment Of Hallux Limitus/Rigidus. *Podiatry Today* Volume 30 - Issue 2 - February (2017)

This content is provided as an educational tool only and is not meant as medical advice in the usage of specific BioPro products. A healthcare professional must use their professional judgment in making any final determinations in product usage and technique. The product's Instructions for Use, should always be reviewed prior to surgery. Postoperative management is patient-specific and dependent on the treating professional's assessment. Individual results will vary and not all patients will experience the same postoperative activity level or outcomes. This information does not constitute medical, legal, or any other type of professional advice and should not be relied upon as such. It is not to be redistributed, duplicated, or disclosed without the express written consent of BioPro, Inc.

Patent No. 8,845,750



BioPro, Inc.
 2929 Lapeer Road, Port Huron, MI 48060, USA
 info@bioproimplants.com | 1-810-982-7777
www.bioproimplants.com

