



Hand & Wrist

Our unique Hand & Wrist portfolio offers proven joint replacements and fixation devices to address common conditions and procedures.

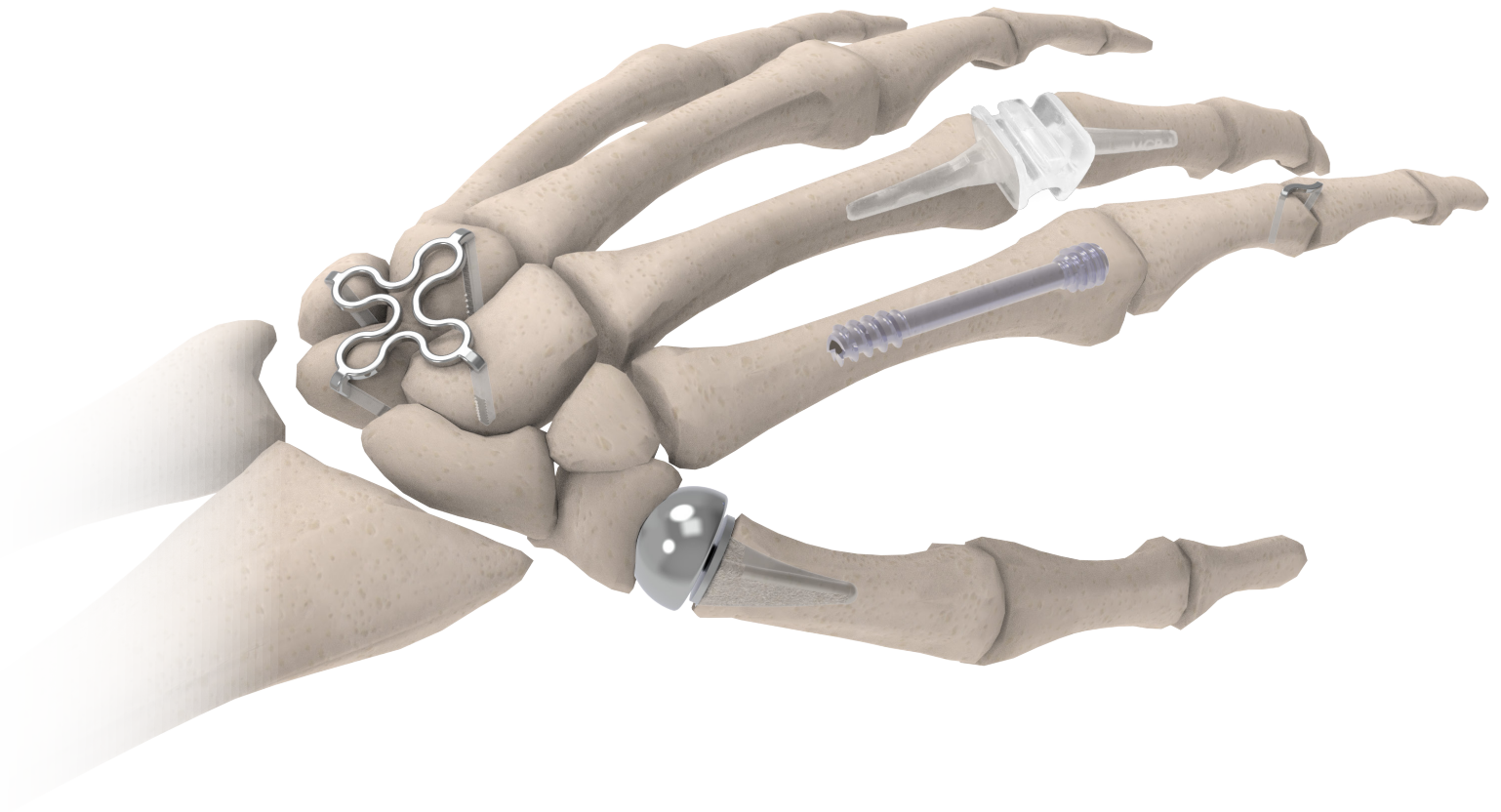

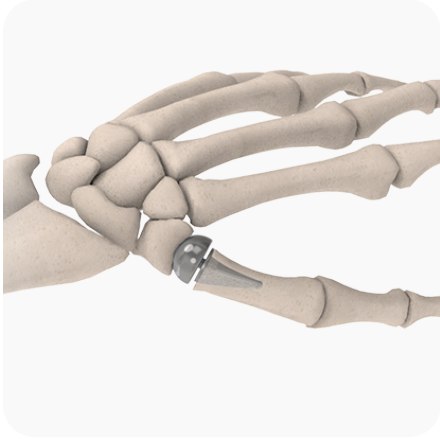

Modular Thumb Implant

A customizable CMC joint implant designed to preserve anatomy, thumb length and pinch and grip strength. The implant is supported by two decades of combined clinical experience and research with published results of 94% survivorship at 6 years.¹



Clover Staple™

3-leg and 4-leg nitinol memory-alloy staple specifically designed for four corner fusions and other limited wrist fusions. The staples activate at patient body temperature, providing dynamic, centralized compression designed to speed up the fusion process.



Digitalis™

The Digitalis Silicone Spacer features a pre-flexed dual hinge design allowing for improved motion and longevity. The spacer is available in 5 sizes for the metacarpophalangeal (MCP) joints and 4 sizes for the proximal interphalangeal (PIP) joints.



Cannulated Compression Screws

A comprehensive line of titanium cannulated screws including headed and headless screws. Go-EZ® headed screws range from 2.0 to 6.5mm while headless designs are available in 2.5 or 3.0mm.



Memory Staple

A 2-leg nitinol memory-alloy staple available in 17 standard sizes. The unique s-bridge staple activates at patient body temperature to provide dynamic compression.

This content is provided as an educational tool only and is not meant as medical advice in the usage of specific BioPro products. A healthcare professional must use their professional judgment in making any final determinations in product usage and technique. The product's Instructions for Use, should always be reviewed prior to surgery. Postoperative management is patient-specific and dependent on the treating professional's assessment. Individual results will vary and not all patients will experience the same postoperative activity level or outcomes. This information does not constitute medical, legal, or any other type of professional advice and should not be relied upon as such. It is not to be redistributed, duplicated, or disclosed without the express written consent of BioPro, Inc.



BioPro, Inc.
2929 Lapeer Road, Port Huron, MI 48060, USA
info@bioproimplants.com | 1-810-982-7777
www.bioproimplants.com

