

Tool VIP Screw System

A screw system specifically designed for minimally invasive bunion surgery.




The System


The TOOL VIP Screw System offers three different diameter screws and specialized instrumentation designed for minimally invasive surgery. The neutral pitch allows surgeons to maintain length of the osteotomy site and optimize stability.




DIAMETER	Ø2.5	Ø3	Ø4
Length	16 to 40 mm	16 to 60 mm	16 to 60 mm
Hexagonal Driver Size	1.7	2.0	2.5
K-wire	Ø1.0	Ø1.0	Ø1.5
Color code	Green	Yellow	Blue

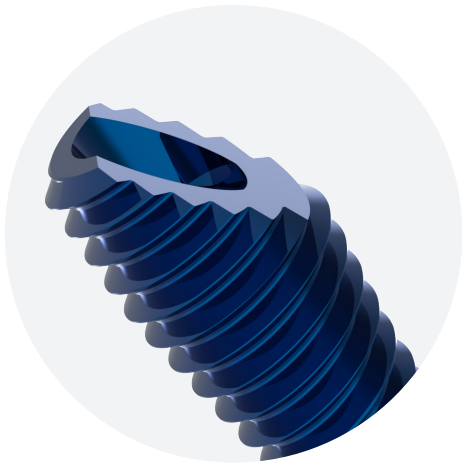

The screws are manufactured from Titanium Alloy and feature a 45 degree beveled head that sits flush with the bone upon insertion.

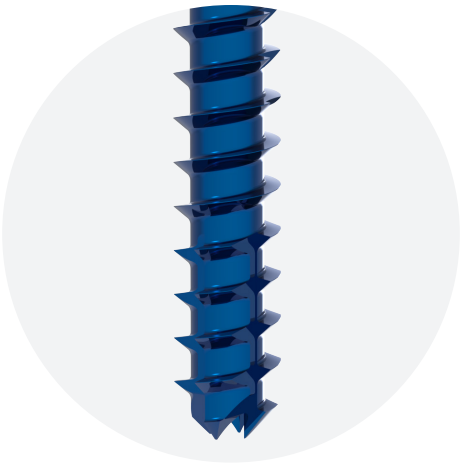

The screws are fully threaded with a neutral pitch and feature self-drilling + self-tapping cutting tips.

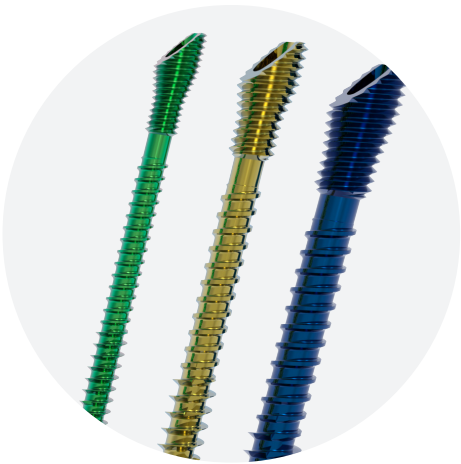

The three different screw diameters allow for use in multiple applications. The screws are color coded and sterile packaged for improved traceability.

Ordering

Length	2.5MM	3.0MM	4.0MM
16MM	SSL125016S	SSL130016S	SSL140016S
18MM	SSL125018S	SSL130018S	SSL140018S
20MM	SSL125020S	SSL130020S	SSL140020S
22MM	SSL125022S	SSL130022S	SSL140022S
24MM	SSL125024S	SSL130024S	SSL140024S
26MM	SSL125026S	SSL130026S	SSL140026S
28MM	SSL125028S	SSL130028S	SSL140028S
30MM	SSL125030S	SSL130030S	SSL140030S
32MM	SSL125032S	SSL130032S	SSL140032S
34MM	SSL125034S	SSL130034S	SSL140034S
36MM	SSL125036S	SSL130036S	SSL140036S
38MM	SSL125038S	SSL130038S	SSL140038S
40MM	SSL125040S	SSL130040S	SSL140040S
42MM		SSL130042S	SSL140042S
44MM		SSL130044S	SSL140044S
46MM		SSL130046S	SSL140046S
48MM		SSL130048S	SSL140048S
50MM		SSL130050S	SSL140050S
52MM		SSL130052S	SSL140052S
54MM		SSL130054S	SSL140054S
56MM		SSL130056S	SSL140056S
58MM		SSL130058S	SSL140058S
60MM		SSL130060S	SSL140060S

Sterile packaged Shannon, Cylindrical and Wedge Burrs are also available for use with or without the Tool VIP Screw System.

Item #	Description
BRB02S62013S	Shannon Std Burr 2.0 X 13mm
BRB01S62008S	Shannon Short Burr 2.0 X 8.0mm
BRB06S62020S	Cylindrical Straight Burr 2.0 X 20mm
BRB07S62514S	Cylindrical Straight Burr 2.5 X 14mm
BRB08S63122S	Cylindrical Straight Burr 3.1 X 22mm
BRB09S63113S	Wedge Burr 3.1 X 13mm
BRB10S64113S	Wedge Burr 4.1 X 13mm





Call us at 1-810-982-7777 to schedule a case today.

The Tool VIP Screw System is a product manufactured by BRM Extremities. For more information and indications for use visit www.brm-extremities.com

This content is provided as an educational tool only and is not meant as medical advice in the usage of specific BioPro products. A healthcare professional must use their professional judgment in making any final determinations in product usage and technique. The product's Instructions for Use, should always be reviewed prior to surgery. Postoperative management is patient-specific and dependent on the treating professional's assessment. Individual results will vary and not all patients will experience the same postoperative activity level or outcomes. This information does not constitute medical, legal, or any other type of professional advice and should not be relied upon as such. It is not to be redistributed, duplicated, or disclosed without the express written consent of BioPro, Inc.



BioPro, Inc.
2929 Lapeer Road, Port Huron, MI 48060, USA
info@bioproimplants.com | 1-810-982-7777
www.bioproimplants.com

